

Series 2800XL


Gymnastic table for therapy acc. to Bobath and Vojta with electric motorized height adjustment

Model 2800-00XLE* white powder coated (RAL 9010)

Model 2800-03XLE* white aluminium powder coated (RAL 9006)

Model 2800-04XLE* grey aluminium powder coated (RAL 9007)

* Exact name of model depending on width of lying surface (see technical data)

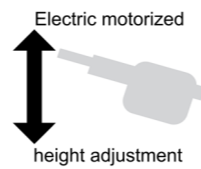

Model 2804-00XLE with options:
Wheel-lifting/fixing mechanism (045N), surrounding foot control (095E), head part stepless adjustable by gas spring (063); upholstery colour: 64-707 chestnut

Basic equipment:



- Lying surface 1-sectioned
- Table not movable, including height adjustment screws to compensate for uneven ground
- Electromotive drive control by foot switch (option: hand switch)
- Integrated safety function for automatic blocking of electromotive drive 
- Freewheeling downwards, automatic declutch in case of resistance
- Dewert-hygiene-equipment 




Table in lowest position



Both pictures.:
Model 2800-00XLE with option:
Movable on twin castors with central lock (046); upholstery colour: 64-717 deep sea

Technical data:

- Electronical height adjustment: 480 mm to 920 mm in 22 sec. by AC-motor with a pressure capacity of 10000 N
- Width of lying surface: **Model 2800XLE:** 1000 mm | **Model 2802XLE:** 1100 mm | **Model 2804XLE:** 1200 mm
- Length of lying surface: 1950 mm
- Load capacity: 250 kg; according to IEC 60601 tested with fourfold static load 

Design:

2 different variants for artificial leather eligible for upholstery covering:

- Skai® Pandoria Plus with velvet-smooth surface
- Skai® Toronto EN with innovative staynu-technology



Available as an option (see pricelist for additional options)

Mobility



Art. no. 040N
Movable on comfort twin castors Ø 100 mm, each with total lock, basic height -15 mm



Art. no. 046
Movable on twin castors, Ø 100 mm with central lock, basic height +20 mm



Art. no. 045N
Wheel-lifting/fixing mechanism with 4 foot levers for controlled lowering

Additional motor control options



Art. no. 090E
Hand switch (low voltage), with magnet on the back to fasten



Art. no. 091E
Foot switch (low voltage)



Art. no. 092E
1 foot switch fixation mounted at underframe



Art. no. 094E
Foot switch rail, lateral, on both sides, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered



Art. no. 095E
Foot switch rail, surround, security adjustment:
Press the switch rail down = the table is lifted
Lift the switch rail up = the table is lowered

Additional safety options



Art. no. B1997E
Foot switch cover, made of brushed VA-steel

Paper roll holder



Art. no. 016K
Paper roll holder, mounted under head part (not protruding), for width 1000 mm

Head part options



Art. no. 063
Head part, upwards and downwards adjustable by gas spring

