



Overview:

An extremely strong and rigid bariatric trolley which can be used as a mobile chair or stretcher for transporting the patient from bed to clinic or to different parts of the hospital. This purpose made trolley has a safe working load and lifting capacity of 320kgs and has electric operation of height and backrest angle to reduce manual handling issues when positioning the patient for transport, various examinations or minor treatment procedures. The dual footrests are manually positioned and have extending foot support boards. The trolley has battery operation for maximum flexibility, and has an integral push handle for steering and transportation.

Benefits


- Lifetime frame & 5 year parts guarantee - see website for further details
- Electric battery operation of chair height and backrest angle for patient positioning
- Flexible rechargeable battery operation complete with charger unit and spare battery
- Electric backrest angle from horizontal to +90°
- The trolley is fitted with fold-down padded side supports with armrests
- The foot sections manually lower from horizontal to -90° and have fold-down foot boards
- Safe working load and lifting capacity of 320kgs (704lb / 50 stone)
- Central locking castors with steering facility to aid transport procedures
- The trolley has an integral transport handle and is supplied with an IV pole
- Low patient transfer height of 53cm (21")
- Wide choice of options, colours and accessories to meet user requirements




Overall Dimensions

Length: 178cm (70") approx
 Width: 80cm (31") approx
 Height: 53cm - 100cm approx
 (21" - 39") approx

Upholstery Dimensions

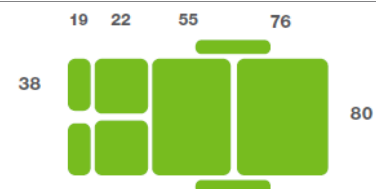

Options & Accessories:

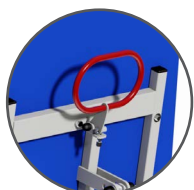

CPR Motor Back



Oxygen Cylinder Holder




D Cushion

NHS Approved Supplier



t: +44 (0)1449 767887
 e: sales@plinthmedical.com
 www.plinthmedical.com

Operational Parameters

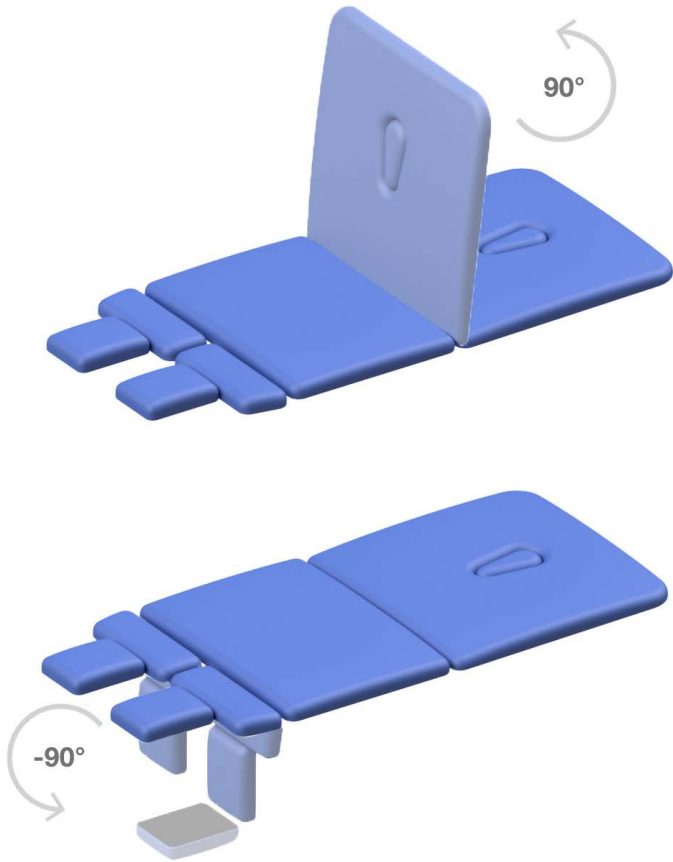

Model BARI 3 - Electric Version

Standard Features

- Designed and manufactured in the UK
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Electric height and backrest operation
- Comfort foam
- Easy operation
- Low running costs
- Hardwearing antimicrobial vinyl
- Safe working load 320kg
- 125mm twin lock castors for patient movement
- Hoist clearance 130mm
- Battery operation complete with charging station and spare battery
- Push handle
- IV pole
- Padded side supports with integral armrests
- Divided legs with fold-down foot boards

Technical Standards

- IP Rating - IPx4
- EEC Classification - MDD93/42/EEC Annex IX Class 1
- Electrical Standard - IEC601 Part 1 + 60601 P (EN)
- Vinyl Flame Retardant to CRIB 5 (BS 7176)
- Foam to BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]
- CE Marked
- Externally audited to ISO9001:2015 by LRQA

Delivery Promise

- Delivery to UK mainland in just 7-10 working days
- Faster service may be available if required
- Overseas delivery 7-10 days plus shipping time

Upholstery Vinyl

- Flame retardant to BS 5852 Ignition Source 5 - CRIB 5 (BS 7176) Medium Hazard
- Independently tested to verify effectiveness at reducing the spread of MRSA, E.coli & other similar bacteria & fungi
- Rated at 500,000 Martindale rubs without damaging or marking the top layer, making it suitable for all high-traffic, heavy-duty areas
- Phthalate free
- Please note: If the upholstery colour is not specified, 'Almond' will be supplied



Jasmine



Battleship



Cobalt



Cocoa



Grape



Mulled Wine



Pillarbox



Gingersnap



Candy



Applemint



Sapphire



Atlantic Blue



Lupin



Denim



Cool Blue



Rainforest



Wasabi



Citrus Green



Marigold



Almond

Plinth Medical Ltd
Barric Lane
Occold, Eye
Suffolk, United Kingdom
IP23 7PX

t: +44 (0)1449767887
e: sales@plinthmedical.com
www.plinthmedical.com

Plinth | Medical